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【M85 Airway Suction Trainer】 11229-060 Oral Cover for training using liquid

Read the instructions below carefully before using.


Setup: Set the Oral Cover at least 5 minutes before starting the training session.
(The Oral Cover just after being set is easier to come off due to lubrication).

Apply the lubricant on the larynx backside of the Oral Cover. (On the circled area of the image above)

Hold the Oral Cover from the front side. Fold it in two by pressing the upper palate with a finger. Open the manikin's mouth carefully and insert the cover into the oral cavity.

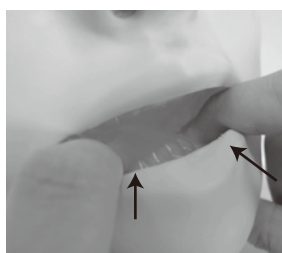Unfold the Oral Cover in the mouth with your fingers.

Hold the Oral Cover's upper and lower ends with your fingers, and carefully remove it.

CAUTION

- ・ Be sure to use the 「lubricant for manikin」 included in the set. Lack or shortage of lubrication may lead to damage the products.
- ・ After drying it well, please store the Oral Cover in a place free of high temperature and humidity.
- ・ Be careful not to have the Oral Cover's surface be in contact with printed objects. The ink might get absorbed by the surface of the Oral Cover, and become impossible to remove.

2016.11.16

■ Main Office and Factory (World Wide)

■ KYOTO KAGAKU CO.,LTD

WEB ● www.kyotogagaku.com

E-MAIL ● rw-kyoto@kyotokagaku.co.jp

15 kitaneokoya-cho Fushimi-ku Kyoto 612-8388, Japan

TEL: +81-75-605-2510

FAX: +81-75-605-2519

■ North and South American regions

■ KKAmerica
Kyoto Kagaku America Inc.

WEB ● www.kkamerica-inc.com

E-MAIL ● info@kkamerica-inc.com

3109 Lomita Boulevard, Torrance, CA 90505-5108, USA

TEL: +1-310-325-8860 (Toll-free in North America: 877-648-8195)

FAX: +1-310-325-8867