KT-27

Elderly Person Simulator Kit

URL ● <http://www.kyotokagaku.com>

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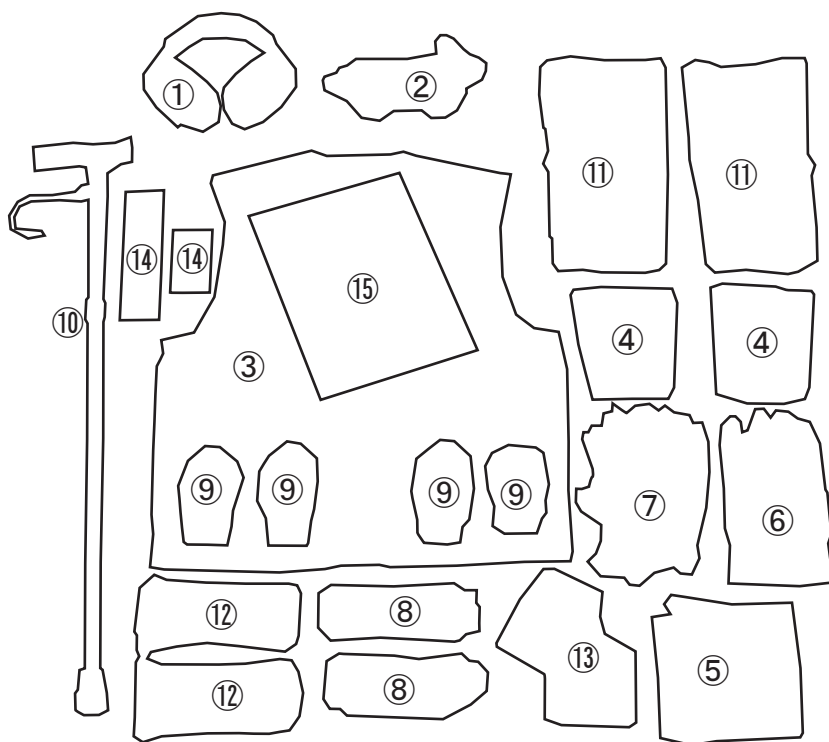
How to wear

Apparatus are to be worn on trainee's dominant foot/leg or hand.


Keep a stable posture when you wear any items.

The product is designed for Educational purpose.


Do not use for any other purpose.

- ① Put the ankle fixture on your dominant foot.
✖Do not take off your shoes when you put on ankle fixture
Use reinforcing velcro belt when the fixture's upper part is too short to wrap the ankle.
- ② Wear the knee fixture on the knee of the same side.
- ③ Put the included weights over the ankle on the same side.
- ④ Wear the vest.
- ⑤ Put the elbow fixture on your dominant hand
- ⑥ Tie the weight around the wrist of the same side.
- ⑦ Wear the goggles.
- ⑧ Wear the earmuff.
- ⑨ Wear the rubber gloves then, wear the white gloves and then hand fixture.
- ⑩ Hold the collapsible stick by other hand and adjust to your size.

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

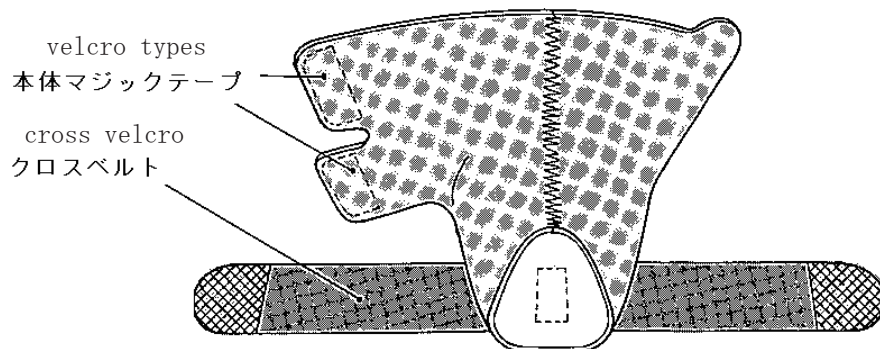
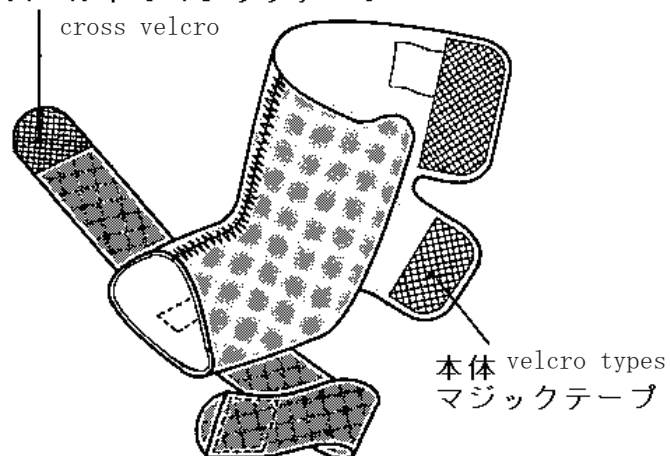
Disengage all velcro tapes for supporter and cross velcro.

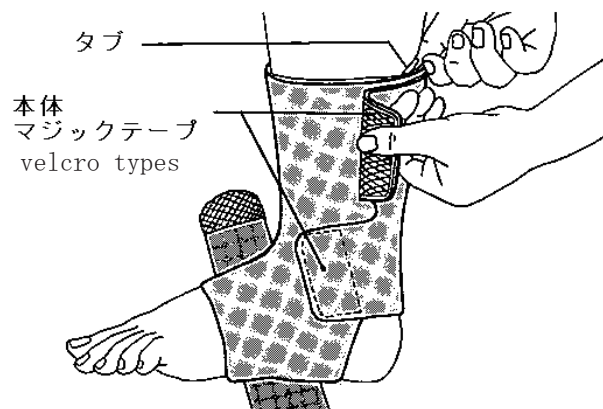
Then, put foot in the hole, wrap the upper part and fix it with two velcro fasteners and cross a velcro on instep.

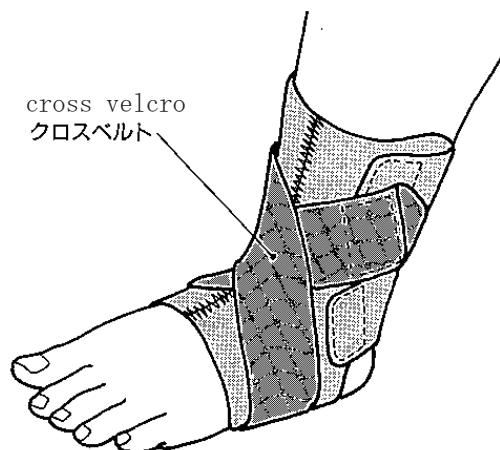
② Knee fixture

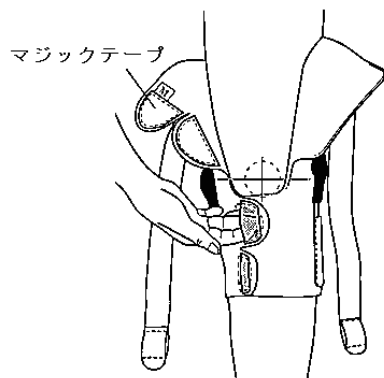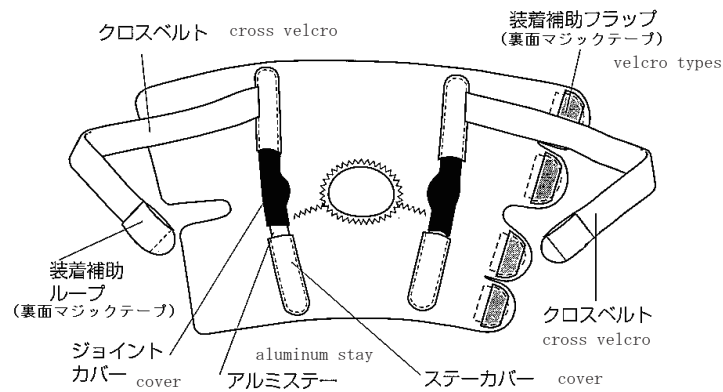
Bend your knee and wear the fixture from the back side of your knee.

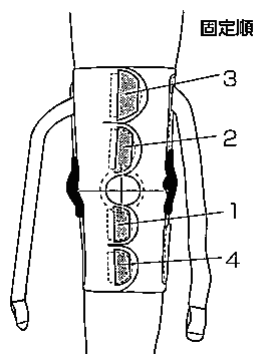Adjust the pad to a middle of your knee, then fasten the velcro fastener in order from 1 to 4.

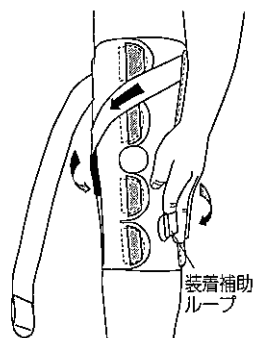※See the illustration.
Cross the velcro belt at the bottom to the top.
Cross the velcro belt at the top to the down.


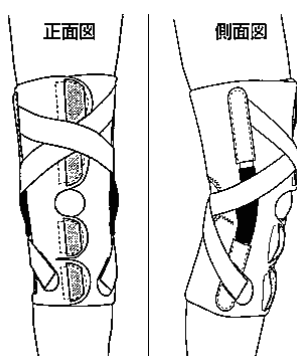
★Caution★

Remove aluminum stay before washing.
Wash with a neutral detergent by hand water temperature should be under 30 degree.
Do not bleache, machine wash, machine dry, iron and dry cleaning.
Dry in the shade.


③ Vest


Put some weights through a Velcro and weight sticks together with type.

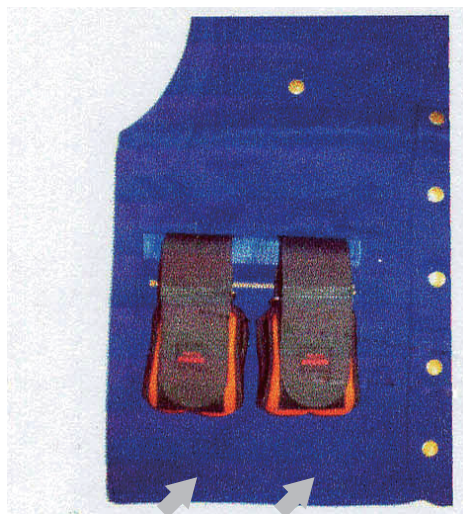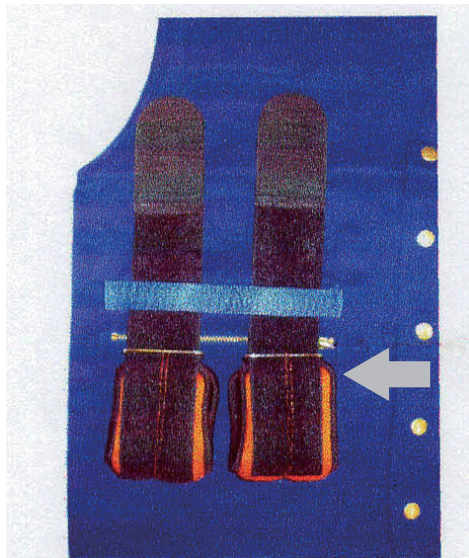④ Elbow fixture

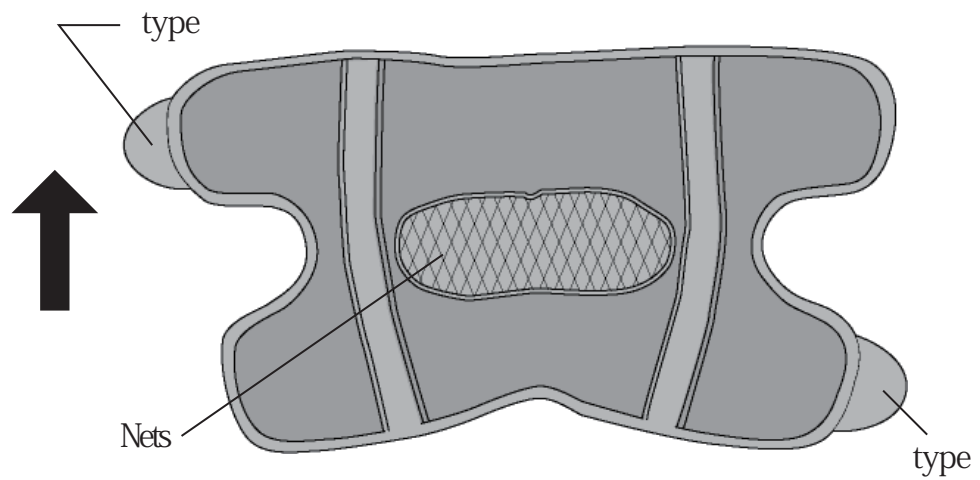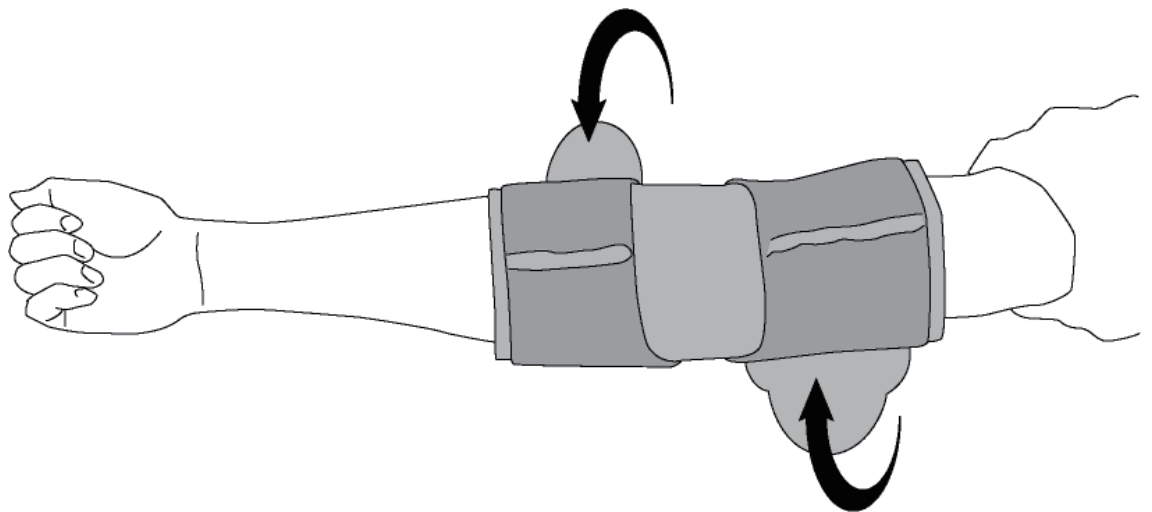Unroll a supporter and put your elbow on it.
Be sure to keep the elbow straight.
Wrap the fixture around the elbow.
Adjust a position of your arm through the
supporter and fasten a type.

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