**⚠ Caution** | Do not mark on the model and other components with pen nor leave printed materials contacted on its surface.  
Ink marks on the models cannot be removed.

# M85 ████████████████████ Airway Suction Trainer

## Contents

- English Manual ..... P.2~P.7

Movie Site

English Site

 **Caution** | Don't mark on the model and other components with pen or leave printed materials contacted on their surface.  
Ink marks on the models will be irremovable.

# M85 Airway Suction Trainer

## Instruction Manual

---

## Contents

---

<b>Introduction</b> .....	P.2
Manufacturer's Note / Dos and Don'ts	
<b>Before You Start</b> .....	P.3
Set Includes	
<b>Before Training</b> .....	P.4
Preparation	
<b>After Training</b> .....	P.5~P.7
Remove the skin on Airway Suction Trainer	
Wash the skin on Airway Suction Trainer	
Reinstall the skin on Airway Suction Trainer	

### Manufacturer's Note

M85 Airway suction trainer is designed to facilitate understanding and training of suction catheter insertion and temporary airway suction.

**Do not use the trainer for any other purpose.**

**Please read the instruction manual before use.**

#### Product Features

1. Three approaches for suction catheter insertion: orally, nasally or through the tracheostomy opening .
2. Temporary airway suction by using simulated sputum.
3. Tracheostomy care and management.
4. Realistic life-size anatomy including mouth, tongue, oral and nasal pharynx, larynx, epiglottis, vocal cords, trachea, and esophagus (facilitates understanding insertion length of catheters).
5. An observation window allows the passage of the catheter to be demonstrated .  
The tips of the provided catheters are colored to help the observation.
6. Trachea and bronchus can be easily flushed clean after use.
7. Face, neck, skin and airway unit can be removed for cleaning and is replaceable.

### DOs and DON'Ts

#### DOs

- Wash the unit after each session.
- Powder the model skin with talcum powder before storage.
- Store the trainer in storage case at room temperature away from heat, moisture, and direct sunlight.

#### DON'Ts

- Do not ever wipe the simulator with thinner or organic solvent.
- Do not pull the skin-airway unit with force. Handle with care.
- Do not mark on the simulator with pen or leave any painted materials in contacted with the surface.  
Ink marks on the model will not be removable.

## Set Includes

Before your first use, ensure that you have all components listed below.


- a. Manikin (life-size, 459\*230\*150mm) . . . . . 1 piece
- b. Catheter (14Fr) . . . . . 2 pcs
- c. Syringe (50ml) . . . . . 1 piece
- d. Tray . . . . . 1 piece
- e. Simulated sputum(150ml×3/1jar) . . . . . 1set
- f. Lubricant . . . . . 1 piece

### Consumables and replacement parts

Product Code	Part Name
11229-070	Simulated sputum/simulated watery stool
11229-020	Replacement skin with trachea for M85 1 piece
11229-040	Catheter for M85 2 piece (14Fr)
11229-050	A bottle of lubricant for M85 1 piece

Simulated sputum

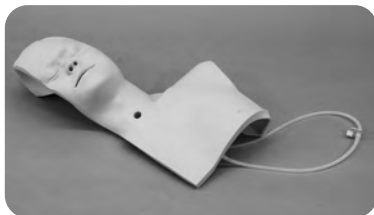Replacement skin with trachea

Catheter

Lubricant

**Please note:** It is recommended to fill the simulated sputum on the provided tray to avoid any stain from the liquid.

### 1 Preparation

1. Open the simulated sputum bag from the notch to open it.



2. Pour the simulated sputum to the accessory bottle.

3. Fill the syringe with simulated sputum.

4. Connect the syringe tip into the three-way tap on the tube from the manikin and inject the simulated sputum into trachea or pharyngeal region.

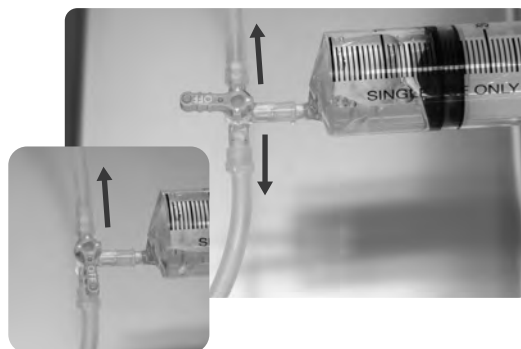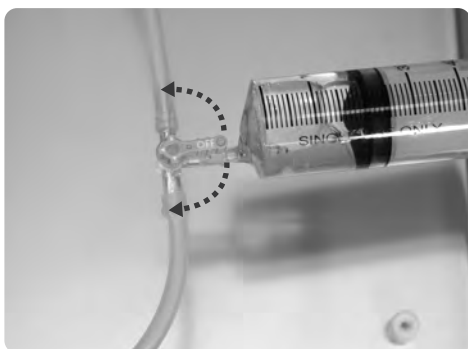5. When the simulated sputum is:


Sticky	Add a little water and stir with stick.
Thin	Stir with stick so as to uniformize the viscosity.

6. Spray the lubricant in the nasal, oral, and tracheostomy openings to facilitate the smooth insertion of catheters.


### 1 Remove the skin on Airway Suction Trainer

1. Close the three-way tap.

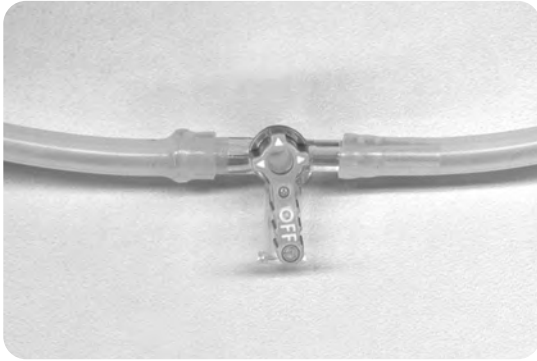2. Detach the side plate by removing the screws.

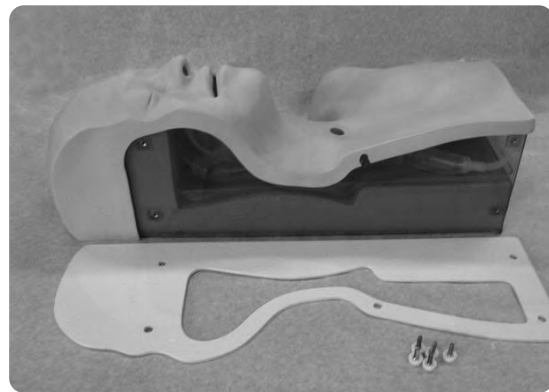3. Remove the skin on Airway Suction Trainer from the base. Be cautious not to spill the liquid from the openings. If the liquid is spilled inside of the manikin, remove the transparent tray from the base and wash it with water.

To remove the skin on Airway Suction Trainer , remove skin from the head side of the model.

### 2 Wash the skin on Airway Suction Trainer

1. Open the three-way tap.
2. Inject the water into the tube by the syringe to eject the simulated sputum.

3. Flush out the inside of the airway unit with running water through each opening.  
Please note: It is recommended to wash the skin after each session.  
This is recommended to avoid drying of the simulated sputum that may create blockage in the tube.

For easy water discharge, detach one of the tube tips from the three way tap.

4. Dry the unit thoroughly. Avoid direct sunlight.

### 3 Reinstall the skin on Airway Suction Trainer

1. Make sure that the nasal tubes are properly placed with the shape of the base.

**Caution**

Do not let ink from pens, newspapers, product manual or other sources contact the manikin. Ink marks on the manikin will be irremovable.

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TEL: +81-75-605-2510

FAX: +81-75-605-2519

**■ North and South American regions:****KKAmerica**

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FAX: 03-3817-8075

The contents of the instruction manual are subject to change without prior notice.

No part of this instruction manual may be reproduced or transmitted in any form without permission from the manufacturer.

Please contact manufacturer for extra copies of this manual which may contain important updates and revisions.

Please contact manufacturer with any discrepancies in this manual or product feedback. Your cooperation is greatly appreciated.

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- For inquiries and service, please contact your distributor or KYOTO KAGAKU CO., LTD.
Kyoto Kagaku America Inc.