

MW67



Caution

Do not mark on the model and other components with pen or leave any printed materials contacted on surface. Ink marks on the model cannot be removed.

Anorectal Manometry Examination Trainer

Instruction manual

Contents



Introduction.....	P.1
Set Includes.....	P.2
Preparation.....	P.3
Training Session.....	P.4
Replacement Part.....	P.5

Introduction

Manufacturer's Note

MW67 Anorectal Manometry Examination Trainer is designed for training in the insertion of a rectal pressure transducer while simulating the anal sphincter in the closed state.

Intended Use

MW67 has been developed exclusively for the training of medical and allied healthcare professionals. Any use other than its intended educational purpose, or any use not in accordance with this manual, is not recommended. KYOTO KAGAKU assumes no responsibility for any accident, injury, or damage resulting from improper use of this product. Handle the product with care. Do not subject it to excessive force, impact, or unnecessary wear. If there are any questions regarding this product, please contact the local distributor or KYOTO KAGAKU. Contact information is provided on the back cover of this manual.

Features

- The **Simulated Pressure Transducer** is printed with scale markings, allowing insertion depth to be checked during training.
- The rectal unit can be removed to visually verify the actual insertion depth of the simulated pressure transducer.


Handling Precautions

• Handle the simulator with care Dropping the simulator or subjecting it to strong impacts may cause damage.	• Do not place printed materials on the product surface. Do not write on the product with a ballpoint pen, permanent marker, or other writing instrument. Ink may be absorbed into the resin surface and become impossible to remove.
• Do not write on the product with ballpoint pens or permanent markers. Ink may be absorbed into the resin surface and become impossible to remove.	• Store away from high temperature, high humidity, and direct sunlight. Improper storage may cause deterioration or damage.
• Discoloration Discoloration may occur during long-term storage or as a result of natural aging. This does not affect product performance.	• Cleaning Wipe the product with water or a mild detergent. If necessary, wipe with alcohol. Never use organic solvents such as thinner.
• Store the product only after it has dried completely. Store in a well-ventilated location after use. Storing the product while wet may result in mold growth.	

Set includes

Set Includes

Verify that all components are present before use.

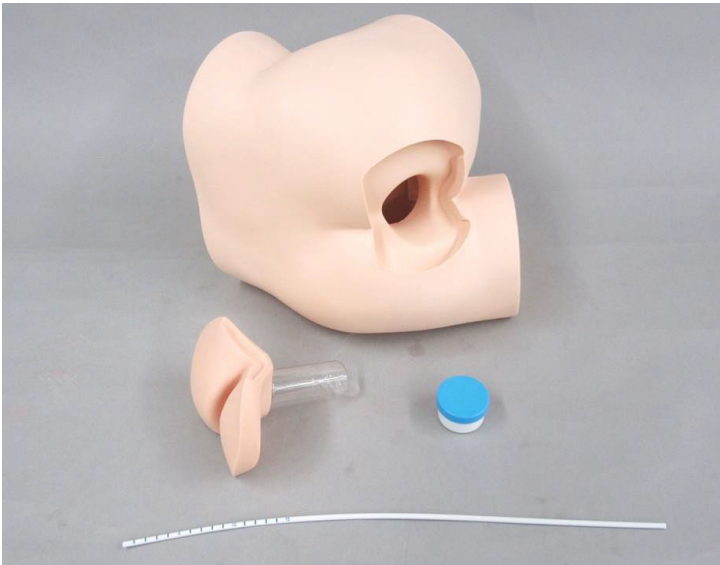

- 1 Lumbar torso
- 1 Rectal unit (pre-installed)
- 1 Simulated Pressure
- 1 Transducer
- 1 Petroleum jelly (Vaseline)


Preparation

Replacing the Rectal Unit


The rectal unit is pre-installed in the lumbar torso.

Replace the rectal unit if the anus or rectum has been damaged.

Insert fingers around the edge of the rectal unit and remove it from the lumbar torso.

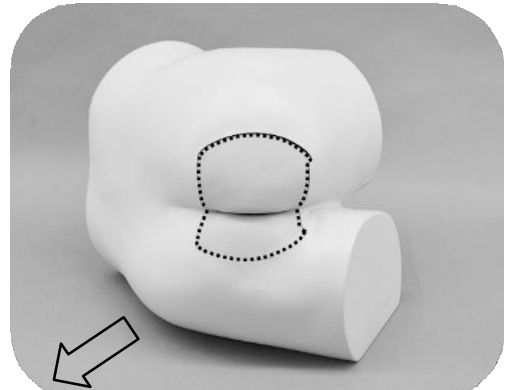

Align the rectal unit with the recess in the lumbar torso, press it into position, and guide the edges by hand until the unit fits flush with the surrounding surface.

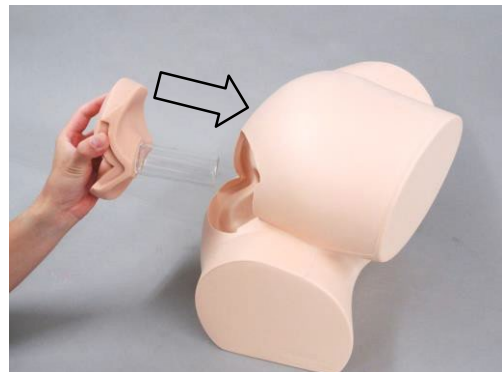


Training session


Insertion of the simulated pressure transducer

Apply petroleum jelly (Vaseline) to the anus before inserting the supplied simulated pressure transducer.

After training, wipe off any remaining petroleum jelly with a well-wrung cloth.

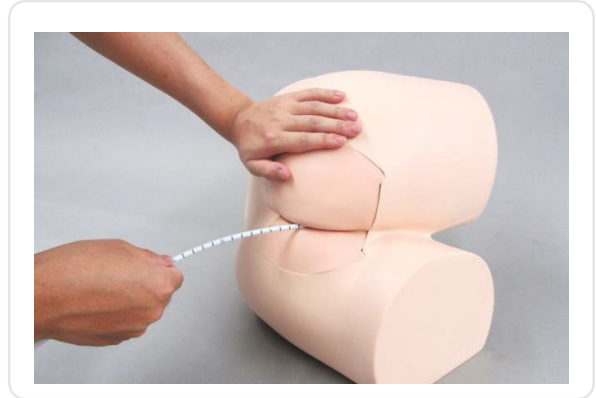

Replacement Part

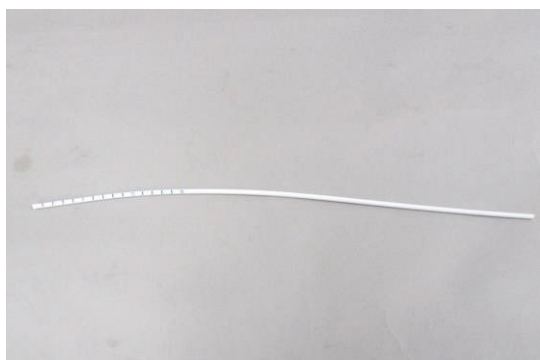

11460-010 Simulated Pressure Transducer

KYOTO KAGAKU co.,LTD

WEB:www.kyotokagaku.com e-mail:rw-kyoto@kyotokagaku.co.jp

Worldwide Inquiries & Ordering

Kyotokagaku Head Office and Factories:

15 Kitanechoya-cho, Fushimi-ku, Kyoto, 612-8388, JAPAN

Tel: +81-75-605-2510 Fax: +81-75-605-2519

All America regions:USA, Canada, and South America



WEB:www.kkamerica-inc.com E-MAIL:info@kkamerica-inc.com

3109 Lomita Boulevard, Torrance, CA 90505-5108, USA

Tel: 1-310-325-8860 Fax: 1-310-325-8867

(Toll free in North America:877-648-8195)

Europe, Russia & Africa

Kyoto Kagaku Europe GmbH.

TEL:+49-69-5060-28160

De-Saint-Exupery-Str.10, 60549 Frankfurt, Germany

The contents of the instruction manual are subject to change without prior notice.

No part of this instruction manual may be reproduced or transmitted in any form without permission from the manufacturer. Please contact manufacturer for extra copies of this manual which may contain important updates and revisions.

Please contact manufacturer with any discrepancies, typos, or mistakes in this manual or product feedback. Your cooperation is greatly appreciated.

July.2027