
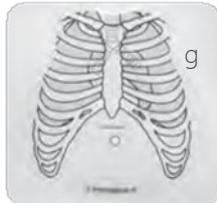
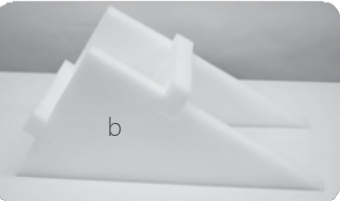


Tube Feeding Simulator Plus


Read the full instruction manual carefully before use.

Set includes

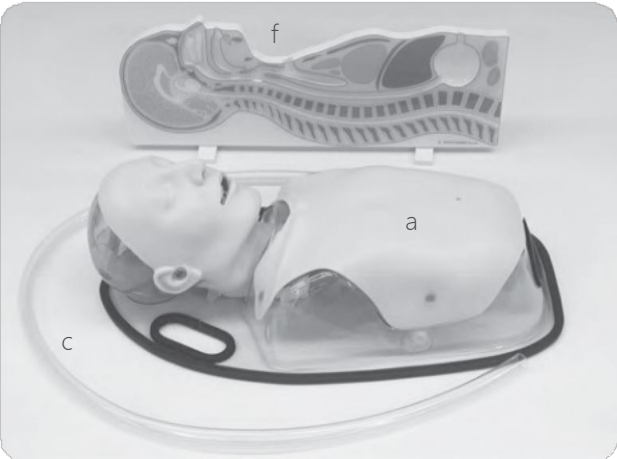


- | | |
|------------------------------|---------|
| a. Male Torso | 1 piece |
| b. Support Base | 1 piece |
| c. Drain Hose | 1 piece |
| d. Funnel | 1 piece |
| e. Plastic cup | 1 piece |
| f. Tube Feeding routes model | 1 piece |
| g. Tube Feeding chest sheet | 1 piece |
| h. Plug for drain port | 1 piece |
| Silicone-based lubricant | 1 |


Preparation (on the bed)

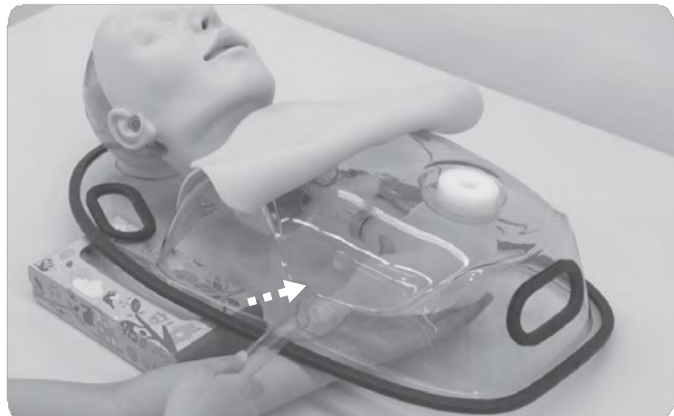

1. Remove the skin on the body and connect the drain hose to the stomach parts.

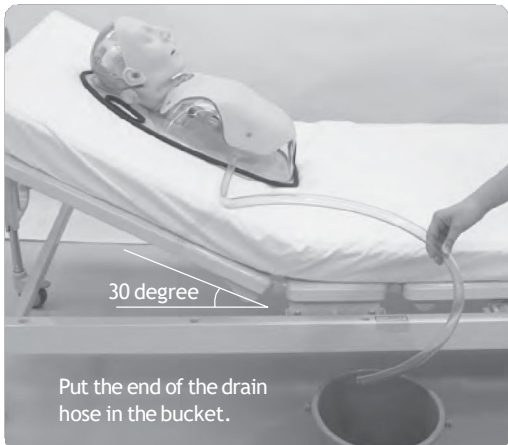

2. Raise up the bed to 30 degree and set the model on the bed. Put the end of the drain hose in the bucket.

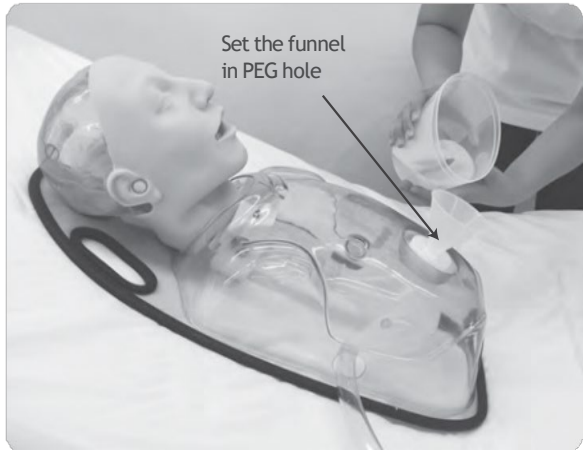


3. Set the funnel in PEG hole and put about 300cc water in the stomach.




4. Set the skin on the body.


Preparation (on the table)

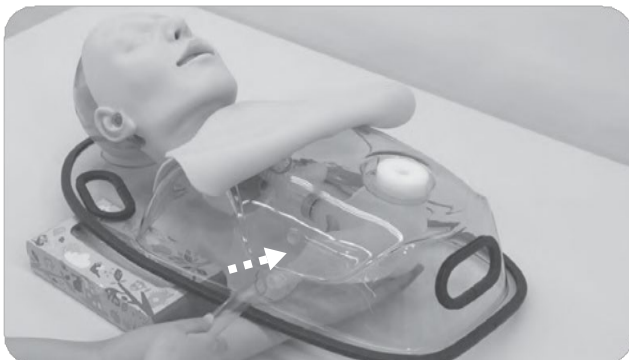


1. Remove the skin on the body and connect the drain hose to the stomach parts.

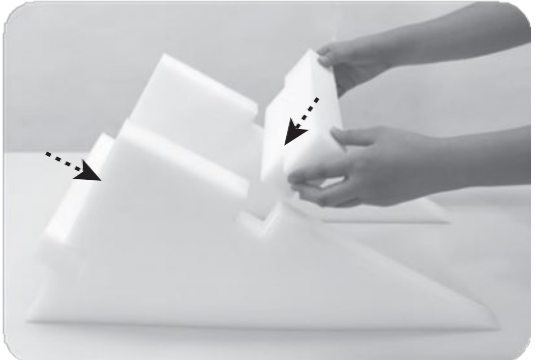

2. Set up the support base.

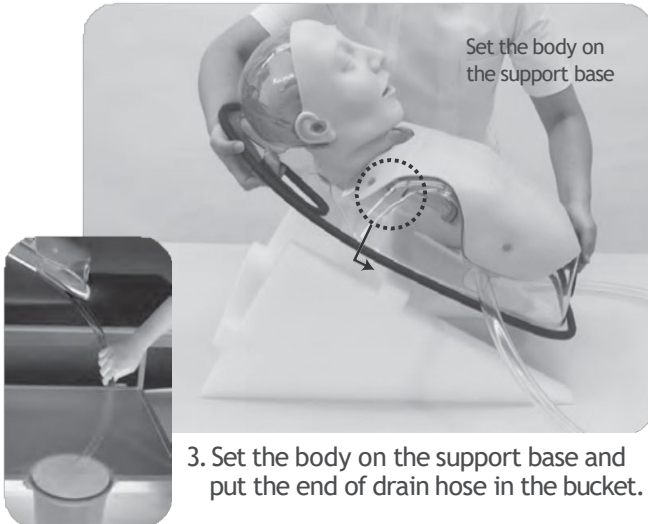

3. Set the body on the support base and put the end of drain hose in the bucket.

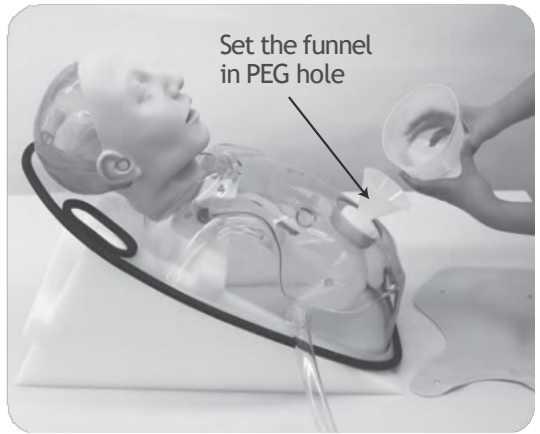

4. Set the funnel in PEG hole and put about 300cc water in the stomach. Then set the skin on the body.

Attention

- Use 14 Fr catheter for this model.
- Apply enough lubricant on the catheter and intranasal passage of the body before starting the practice. (Use attached lubricant)
- Fix PEG tube with about 5ml air. (In case of NIPRO 20Fr catheter)
- Each feeding of a package nutrition, drain the content in the stomach. And add 300ml water again for next feeding.
- Do not leave tapes on the model for long time.

KYOTO KAGAKU co.,LTD

■ Main Office and Factory
15 Kitanechoya-cho Fushimi-ku Kyoto 612-8388, Japan
Telephone : 81-75-605-2510 Facsimile : 81-75-605-2519

URL ● <http://www.kyotokagaku.com>
e-mail ● rw-kyoto@kyotokagaku.co.jp

■ LA Office (for USA , CANADA and Mexico customers)
3109 Lomita Boulevard, Torrance, CA 90505 ,USA
Telephone : 1-310-325-8860 Facsimile : 1-310-325-8867