


PH-7 4 Bone Scintigraphy Quality Assurance Phantom

ISMMType (SIM bone Phantom)

Instruction Manual


Contents

Before you start	
Product contents and precautions.....	P.1
Specifications.....	P.2
Preparation	
How to set up.....	P.3

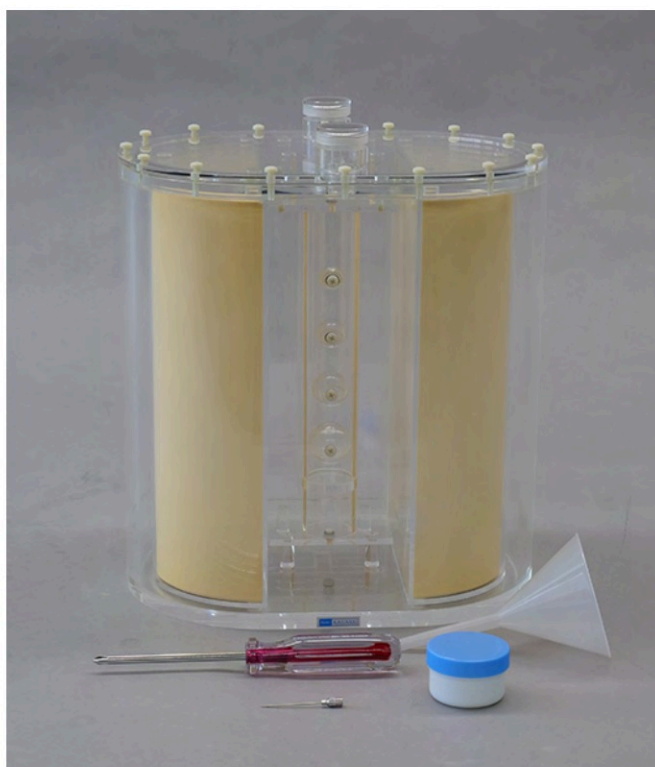

Set Includes

Before your first use, ensure that your product includes all the components listed below.

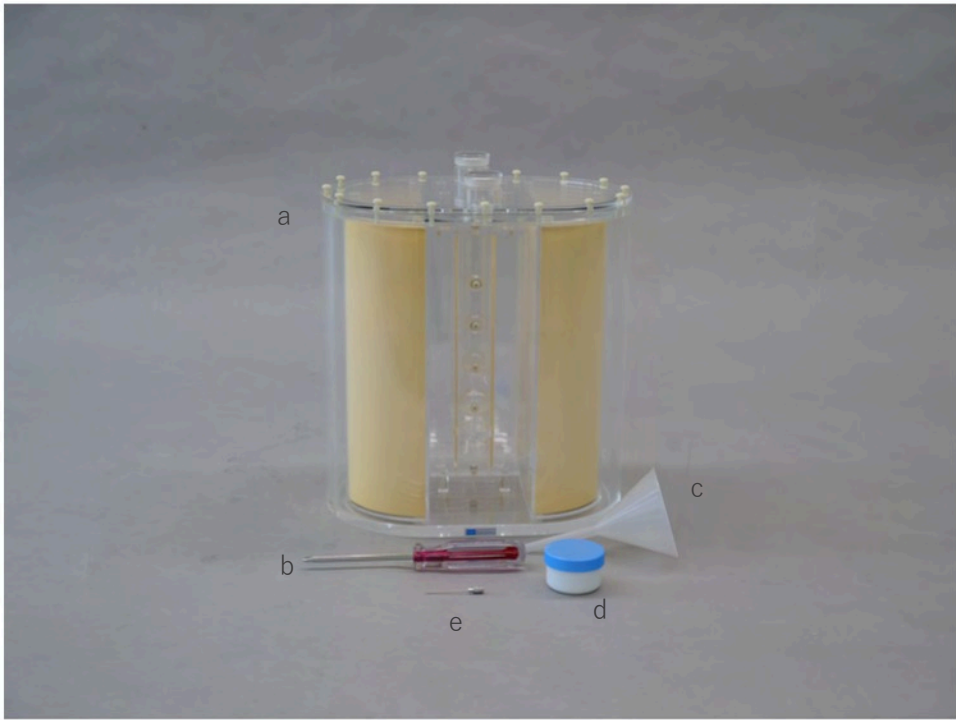

a Phantom body (Including left and right lung fields and bones)	1	d Petroleum jelly	1
b Screw driver	1	e Needle	1
c Funnel	1		

 Precaution

Please handle with care. The phantom
is made of hard resinous material; the phantom can be damaged if dropped or hit by a hard object.

Please clean the phantom with water or
pH neutral detergent. Please do not use organic detergent such as thinner.

Please do not store in high
temperatures, high humidity or in direct sunlight. It may lead to deformation or damages.

Please do not write with a pen.
Writing directly on the phantom with felt-tip pen or permanent marker will not wash off.




Phantom's specs (External dimensions, material)

W310 X D210 X H355 mm (Including protrusions)


Outer container/vertebral body:···Acrylic resin

Inner lung field: ··· tough lung (PVA acetal product)

Drawing

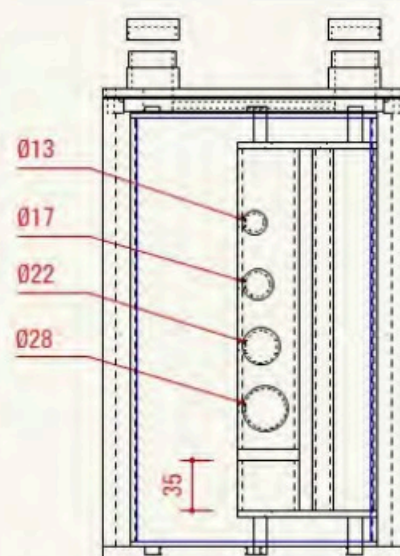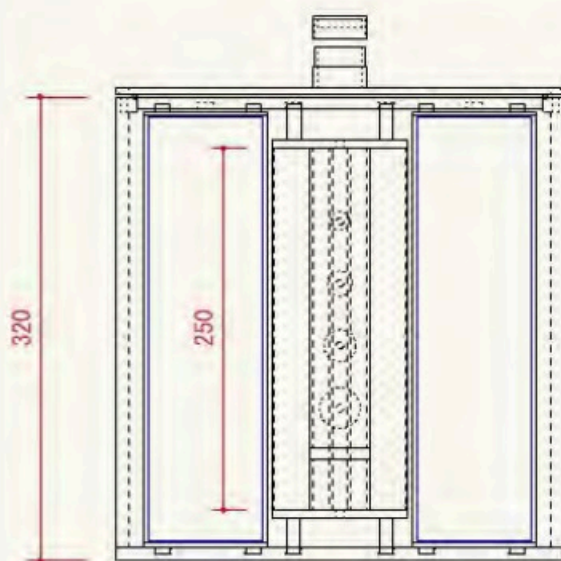
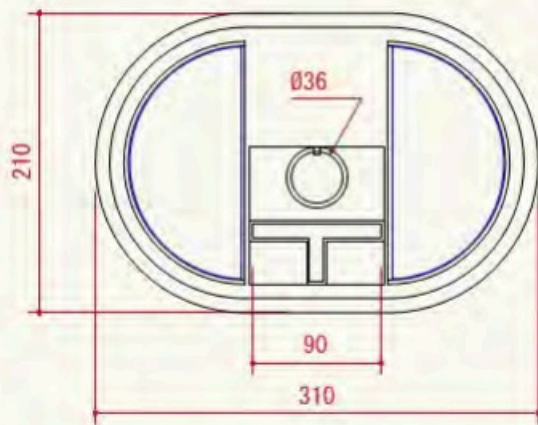


Detachment of the lung field

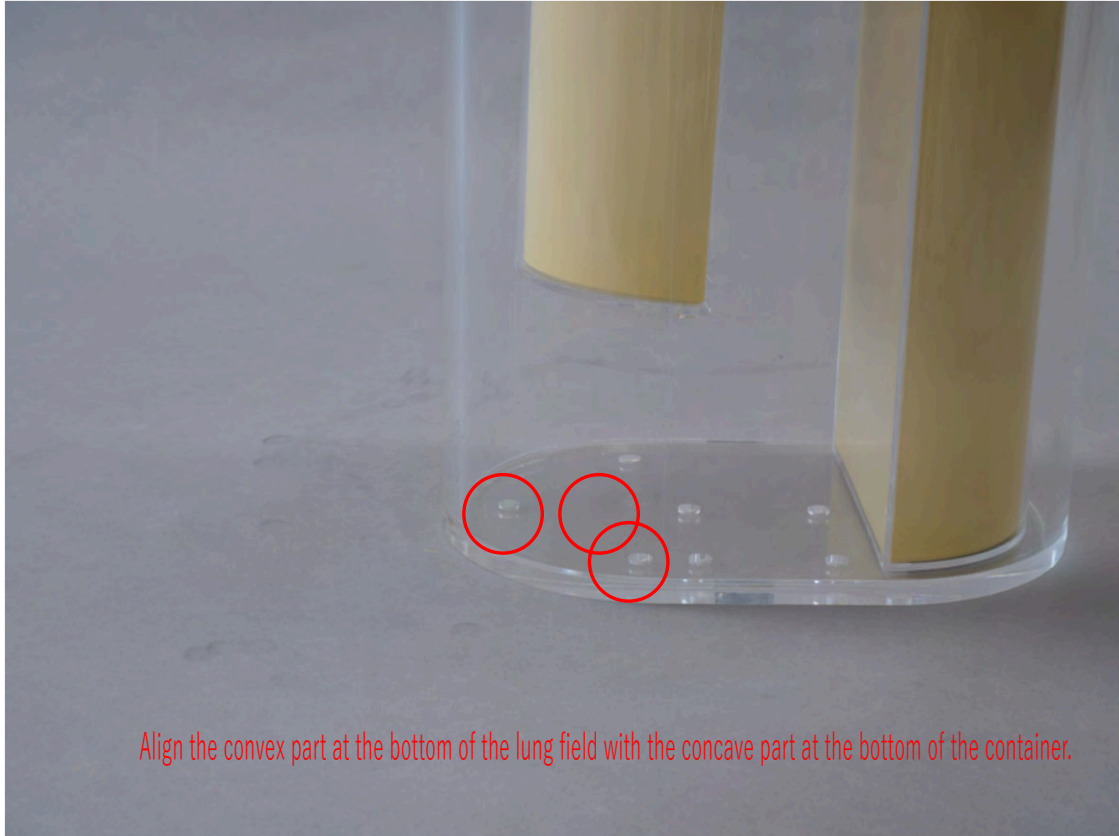

*The container can be used as a chest with the lung fields attached, or as an abdomen without the lung fields.

Assembling the spine

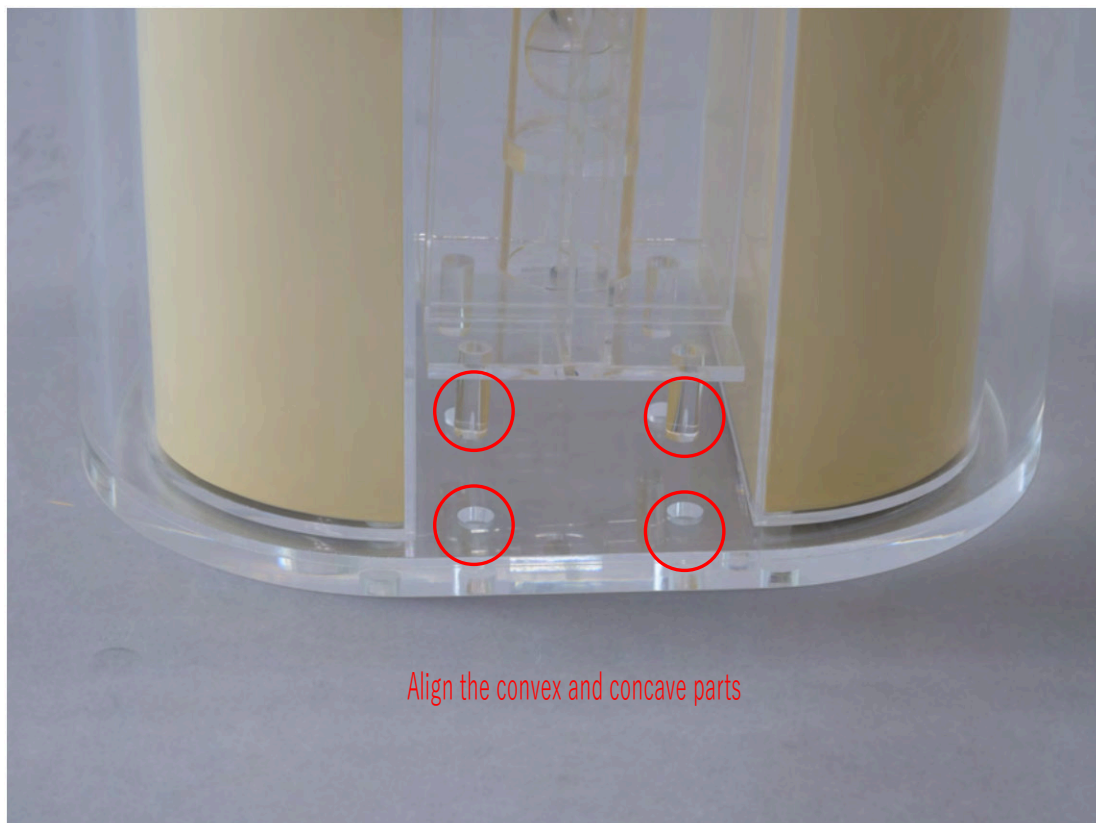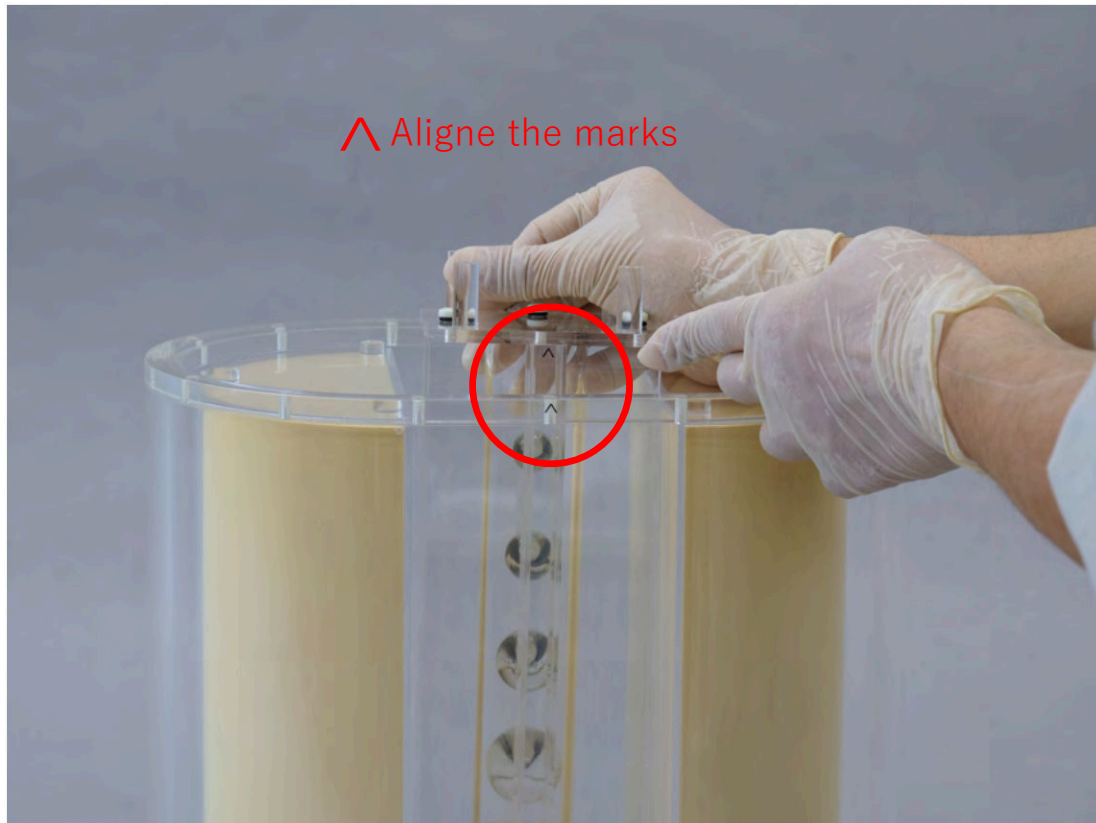


Fixing the lid



Sim² Bone Phantom Volume measurement record

			Lot No.	KYOTO KAGAKU CO.,LTD

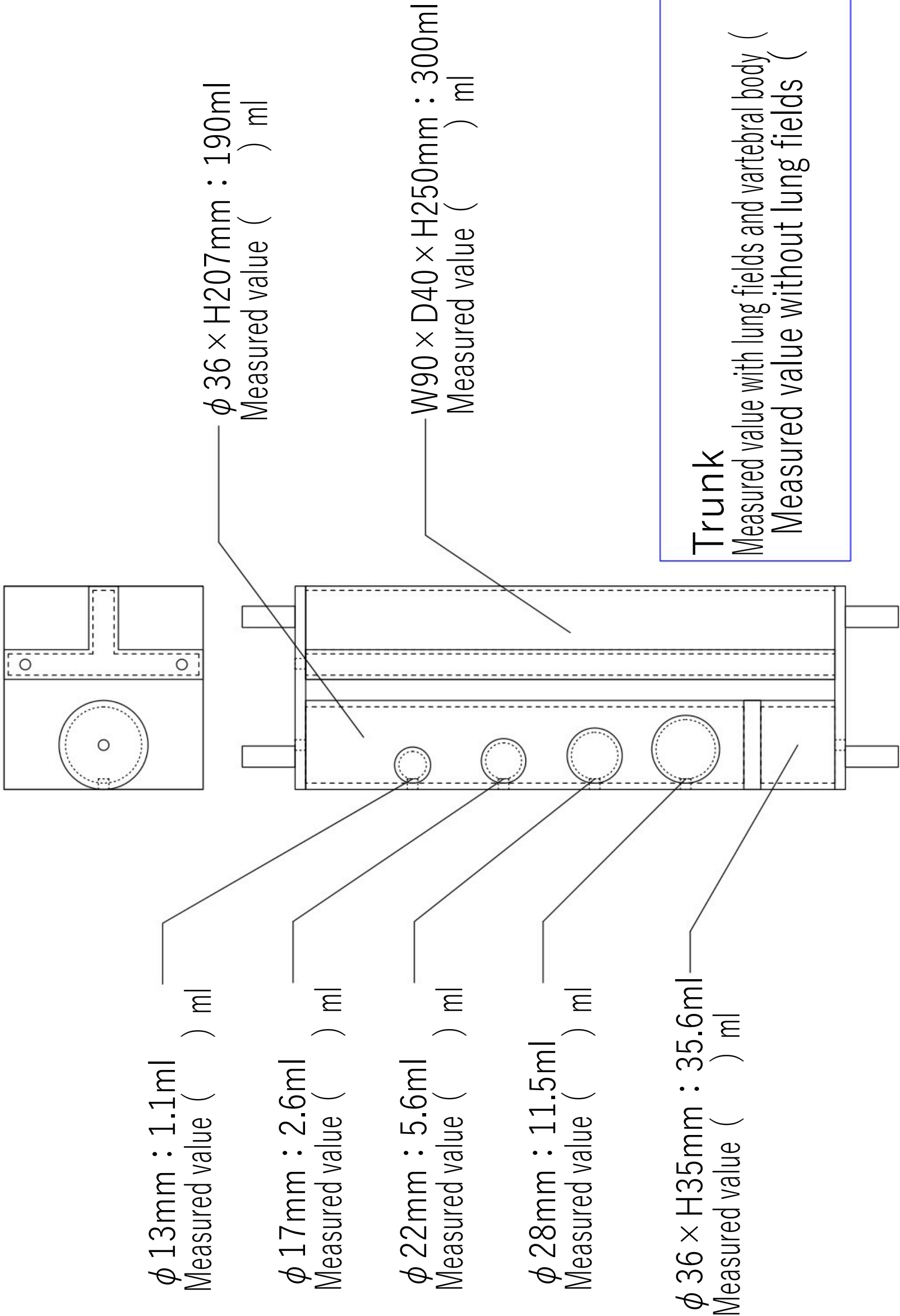

KYOTO KAGAKU CO.,LTD
15 Kitanechoya-cho, Fushimi-ku, Kyoto, 612-8388, Japan

URL <http://www.kyotokagaku.com>
e-mail rw-kyoto@kyotokagaku.co.jp