


MW7B

NEW

Blood Pressure Measurement Trainer ver.2

Become familiar with the physiology of hypertension by blood pressure observations


With 9 preset cases, the simulator allows student lead practice in accordance with the blood pressure data obtained

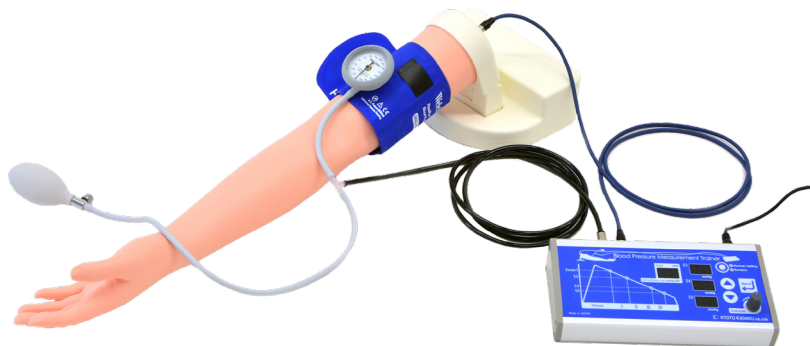

Manual Inflation



Korotkoff Sounds

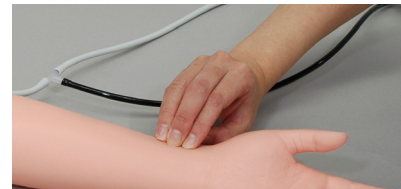

Palpation

FEATURES

- 1 | Nine preset cases:
- 2 | Blood pressure values can be digitally controlled at three points - Phases I, IV and V
- 3 | Realistic auscultation of phase I to V Korotkoff sounds
- 4 | Radial artery pulses synchronized with sound and cuff pressure
- 5 | Deflation rate can be numerically monitored
- 6 | Excessive cuff pressure is warned by beeping
- 7 | Random play mode for self-learning

SKILLS

- | Cuff placement
- | Manual cuff inflation
- | Auscultation of Korotkoff sounds
- | Palpation of radial pulses
- | BP reading
- | Cuff deflation

👍 Nine Preset Cases

- 1 | hypotensive
- 2 | normal
- 3 | prehypertension
- 4 | mild hypertension
- 5 | moderate hypertension
- 6 | severe hypertension
- 7 | isolated systolic hypertension
- 8 | auscultatory gap
- 9 | absence of Phase V Korotkoff sounds

DESCRIPTIONS

SPECIFICATIONS

Size: 67x18x22cm/26.4x7.1x8.7 inch Weight: 2.5kg/ 5.5lbs

Power: AC100V-240V 50/60Hz

Power Consumption: 10W

MATERIALS

Soft resin, Latex free

SET INCLUDES

| | |
|-----------------|-----------------------------------|
| 1 arm model | 1 aneroid sphygmomanometer unit |
| 1 controller | 1 instruction manual |
| 1 AC adaptor | 1 carrying bag (no photo) |
| 1 power cable | |