

MW49

# Simulator Intravenous ARM III

**M50B-B has been renovated to meet a wider range of training requirements.**



## FEATURES

- 1 | Two types of injection pads for venous line placement (forearm/ back of hand)
- 2 | Two types of injection pads for cubital fossa (regular and durable)
- 3 | Arm vessels anatomy sleevelets for visual understanding

## SKILLS

- | Anatomy of the arm
- | Vacuum tube blood collection
- | Palpation /Intravenous injection
- | IV cannula placement and drip infusion
- | Blood collection with a syringe/Flash- back confirmation

### Two types of injection pads for cubital fossa - works with circulation pump




Pad for back of hand



Pad for ventral / dorsal of forearm



### Four areas of puncture sites

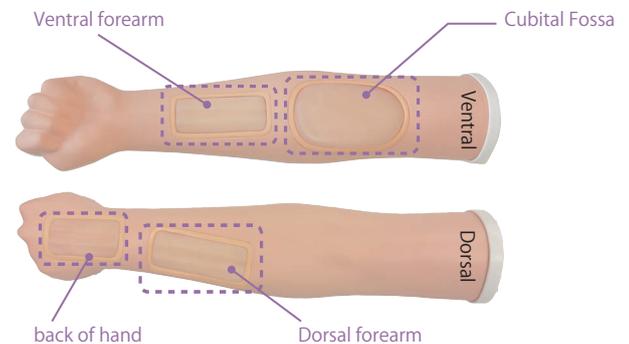

### Two types of injection pads for venous line placement - works with bottles/stands



#### Regular pad

Realistic sensation. Texture, sensation and resistance to the needle. Three accessible venous lines represent basilic vein, median vein, and cephalic vein.

#### Durable pad

Durable training pad to maximize training opportunities. Long-life tissue materials. Both faces of the pad enclose two venous lines for training.

## DESCRIPTIONS

### SPECIFICATIONS

Size : W25 x D11x H3 cm each Power : AC100-240V 50/60Hz  
Power Consumption : 5W

### SET INCLUDES

2 arm model	1 circulation pump with a plastic jar
2 supporting stand	2 bottle with stand (for sim blood)
2 injection pad (R pad)	10 cotton swabs red
2 injection pad (D pad)	2 arm vessels anatomy sleevelets
4 injection pad (forearm/antebrachial)	1 storage case
2 injection pad (back of hand)	1 instruction manual

### REPLACEMENT PARTS

11267-300-03 R pad (set of 2)
11430-020 D pad (set of 2)
11388-200 pad for forearm/antebrachial (set of 2)
11388-300 pad for back of hand (set of 2)
11267-300-02 supporting stand
11267-300-05 circulation pump
11388-400 10 cotton swabs red
11430-010 arm vessels anatomy sleevelet

### RECOMMENDED DEVICES

- IV cannula 23G
- syringe needle 23G
- butterfly needle 23G