

M74



Gluteus Maximus Injection Simulator

IM training for the three ways of injection

This simulator is designed to practice for the three ways of injection into the hip and the buttocks. The incorporated indicator shows the accuracy of needle placement and its depth.

Product Supervision

Dr. Hiroshi Koyama,
Medical Director, General Internal Medicine
National Hospital Organization Kyoto Medical
Center



FEATURES


- 1 | The life-sized lower torso manikin consist of three layers: skin, muscular and skeleton.
- 2 | Greater trochanter, anterior superior iliac spine and scrum can be palpated to identify the injection sites
- 3 | Three methods for intramuscular injection: QuarterTrichotomy, Hochstetter Method and Clark Method.
- 4 | The indicator shows the accuracy of the needle placement and its depth.
- 5 | Correct needle placement can be confirmed by being able to inject liquid.
- 6 | No needle trace remains.

SKILLS

| Intramuscular injection



Hochstetter Method

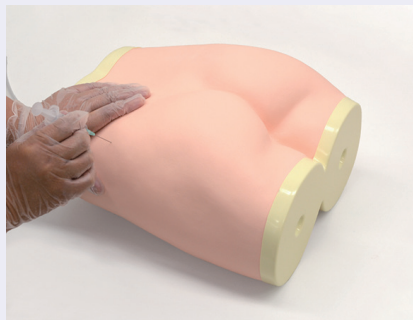

Clark / Sachtleben Method



DESCRIPTIONS

SPECIFICATIONS

Size: W32 x D35 x H20 cm Weight: 3.2kg / 7.05lbs
W20 x D16 x H13 cm

SET INCLUDES

1	lower torso manikin	1	instruction manual
1	stand for Hochstetter method		
1	indicator		
1	power AC adapter		
1	set of sensor pads		

MATERIALS

Soft resin
Latex free

REPLACEMENT PARTS

11213-010	skin
11213-020	sensor for puncture pad
11213-030	sensor for Hochstetter's site
11213-040	sensor for Clark's site

RECOMMENDED DEVICES

IM needle 21-23G