


Cervical Assessment Simulator

Realistic yet durable material for repetitive practice of cervix scoring



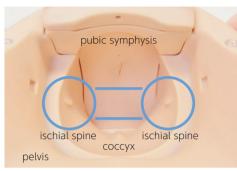

Station zero, the line that connects ischial spines can be recognized.

FEATURES

- 1 | Genital area made of innovative material that excels in elasticity, durability and resilience offers repetitive training.
- 2 | 5 types of cervical dilation models and fetal descent positions can be interchanged with a single action.
- 3 | Anterior/posterior fontanelles to determine the fetal position

SKILLS

- | Cervix scoring in the first stages of labor
- | Check fetal head position
- | Determination of bishop score (cervix score)
- Speculum examination (i.e. Pap smear testing)
- | Urinary catheterization
- | Cervical cytology

KEY FEATURES



Cervical dilation insert can be attached by a single action.



Descent of the fetus can be set by the lever on the holder (from -3 to +3).

DESCRIPTIONS



External genitalia that excels in elasticity, durability and resilience

SPECIFICATIONS

Size: W74 x D44 x H33.5 cm W29 x D17 x H13 inch

Weight: 10.3 kg / 22.7 lbs

SET INCLUDES

- 1 maternal body
- genitalia unit
- vagina unit
- set of cervical dilation inserts (5 variations)
- 1 holder for the cervical dilation inserts
- - 1 lubricant (500 ml)
 - talcum powder instruction manual
 - towel

MATERIALS

Soft resin, hard plastic Latex free

REPLACEMENT PARTS

11415-010 1 set of cervical dilation inserts (5 variations)

11415-040 1 lubricant (500 ml) 11415-030 1 genitalia unit

(for pelvic examination and delivery assistance)

RECOMMENDED DEVICES

OPTIONAL PARTS

11416-020 placenta

11416-100 vaginal delivery module 11418-100 uterine involution module