




# Laparoscopic Suture Training System -Laparo Cube

**Quantitative Feedback combined  
with Realistic Hands-on Experience**

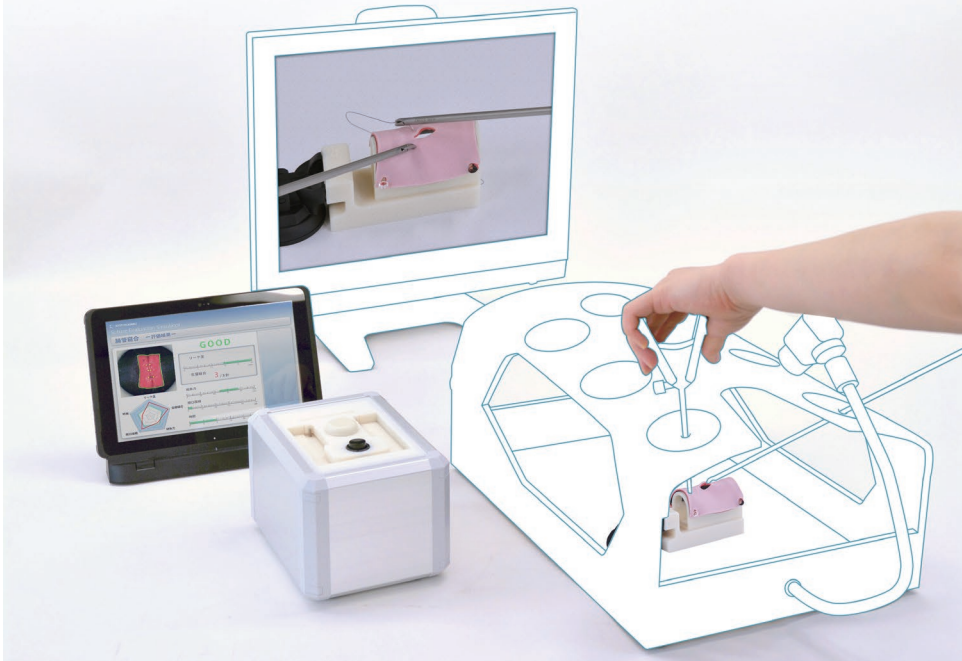

## FEATURES

- 1 | Innovative assessment criteria and quantitative measurement provide objective feedback to trainers.
- 2 | Trainees' personal skills assessment data can be saved for reviewing and analysis.
- 3 | Actual clinical instruments can be used.
- 4 | Realistic and lifelike suturing sheet.
- 5 | Able to create your own passing standard or adjust criteria threshold to trainees' level.


## SKILLS

- | Laparoscopic suturing (simple interrupted suture)

## KEY FEATURES

### Judging success or failure

- 1 | Is the wound sutured without a leak?
- 2 | Is suture completed through muscular layer to mucosal layer?
- 3 | Is ligature appropriate for suture?
- 4 | How much open area is there in mucosal layer?
- 5 | How long does it take for suturing procedure?

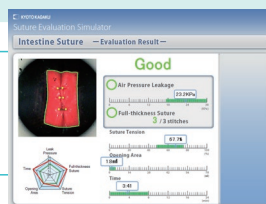

### Life-like Intestinal Tract

Intestinal tract has four layers with respective characteristics. When holding with forceps and suturing, these four layers react like real tissue.

Easy to set up!



## DESCRIPTIONS

### SPECIFICATIONS

Size: W15 x D20 x H13 cm /  
W5.9 x D7.9 x H5.1 inch

Weight: 1.5kg / 3.3lbs  
(Assessment unit)

Power: AC100V-240V  
Power Consumption:  
26VA : Assessment unit  
40VA : PC

50/60Hz

### MATERIALS

Soft resin  
Latex free

### RECOMMENDED DEVICES

-

### SET INCLUDES

1 laparoscopic suturing assessment unit	1 PC
1 power cable for assessment unit	1 power cable for PC
1 set of skin sheet holder	30 pieces of suturing sheet (intestinal tract)
1 cable	1 container
1 USB Cable	1 instruction manual

### REPLACEMENT PARTS

11395-010 30 pieces of suturing sheet  
※1 set consist of inner layer& outer layer

Box trainer, display monitor, and forceps and camera are not included in set