


Blood Pressure Measurement Trainer II

Become familiar with the physiology of hypertension through blood pressure observations

With 9 preset cases, the simulator allows student-led practice in accordance with the obtained blood pressure data

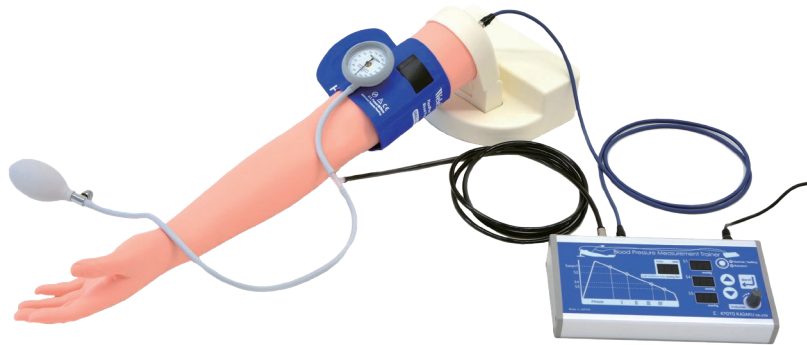

Manual Inflation



Korotkoff Sounds



Palpation

FEATURES

- 1 | Nine preset cases
- 2 | Blood pressure values can be digitally controlled at three points
- Phases I, IV and V
- 3 | Realistic auscultation of Korotkoff sounds phase I to V
- 4 | Radial artery pulses synchronized with sound and cuff pressure
- 5 | Deflation rate can be numerically monitored
- 6 | Excessive cuff pressure is warned by beeping
- 7 | Random play mode for self-learning

SKILLS

- | Cuff placement
- | Manual cuff inflation
- | Auscultation of Korotkoff sounds
- | Palpation of radial pulses
- | BP reading
- | Cuff deflation

Nine Preset Cases

- 1 | hypotensive
- 2 | normal
- 3 | prehypertension
- 4 | mild hypertension
- 5 | moderate hypertension
- 6 | severe hypertension
- 7 | isolated systolic hypertension
- 8 | auscultatory gap
- 9 | absence of Phase V Korotkoff sounds

DESCRIPTIONS

SPECIFICATIONS

Size: W70 x D19 x H17 cm / W27.7 x D7.5 x H6.7 in. Weight: 2.65kg/5.84lbs
Power: AC100V-240V 50/60Hz
Power Consumption: 10W

MATERIALS

Soft resin, Latex free

SET INCLUDES

1 arm model	1 aneroid sphygmomanometer unit
1 controller	1 instruction manual
1 AC adapter	1 carrying bag
1 power cable	1 quick start guide
1 connection cable	